


## Rocking Table KICK II

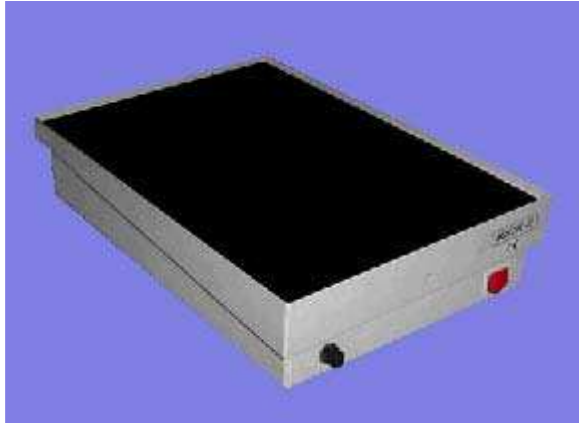

### Rocking Table KICK II:

Thermostatable large size incubation platform, suitable for various incubation and preparation techniques (e. g. ELISAs, immunodiffusion, immunoblotting, histochemistry, in-situ-hybridization), cell- and microbiological culture experiments and for the accelerated processing of staining/destaining procedures with electrophoretic gels and blotting membranes.

- A rugged and robust heavy-duty lab instrument with loading capacity up to 5 kg and an effective area of 45 x 30 cm covered with an anti-slip support.
- Constructed for continuous operation.
- Electronic regulation of lifting frequency (10 - 100/min) with newly developed lifting technique (lifting angle 5°, lifting height 3 cm) performing a continuous light "kick" agitation to prevent gels and membranes from sticking to their containment support!
- Optional doubling of capacity by additional platform attachment (45 x 30 cm) with anti-slip support
- Optional heating and cooling by attachment of up to two HOMOTHERM thermoplates fitting into the resp. platform.

### Ordering information

| Product               | Format     | Cat. No. |
|-----------------------|------------|----------|
| Rocking Table Kick II | 45 x 30 cm | KICK II  |