


## OMX-S

### Peptides Ready for MS in only 1 Hour!

The new OMX-S tube is a multifunctional device developed by OMX, a proteomics company located in Munich/Germany. The sophisticated design allows to carry out all steps, picking, digesting, sample preparation for MS in the same tube. The benefits are no loss of sample material, easy to handle, no cross-contamination.

And most appealing - savings in time of up to 90 % are routinely obtained compared to the standard EMBL protocol.

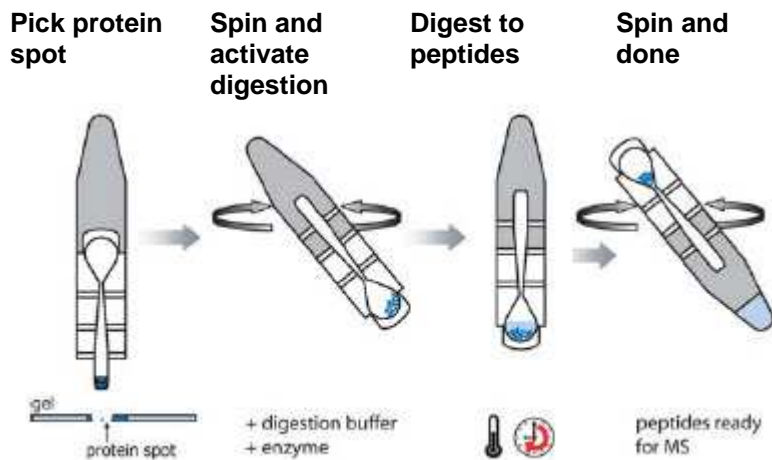

- Multifunctional device, easy to handle
- Fast, savings in time up to 90 % (referring to the EMBL standard protocol)
- No loss of sample material
- No cross-contamination

### OMX-S: Multifunctional device, easy to handle

#### Specifications

##### Material

Polypropylen and glass

##### Diameter of picker


1.8 mm

##### Reactor Volume

max. 40 µl

##### OMX-S length

50 mm

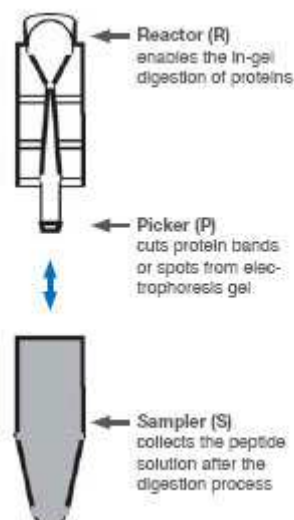

#### Ordering Information

Product	Quantity	Cat. No.
OMX-S	12 pieces	39240.01
Device	2 pieces	39240.02