

SERVA Native Marker Liquid Mix for BN/CN

Kat.-Nr. 39219


Product Description:

General Information

Ready-to-use protein standard for Blue Native (BN) and Clear Native (CN) gel electrophoresis containing 6 different proteins with a molecular weight range between 21 and 720 kDa.

Application

- Blue Native and Clear Native gel electrophoresis

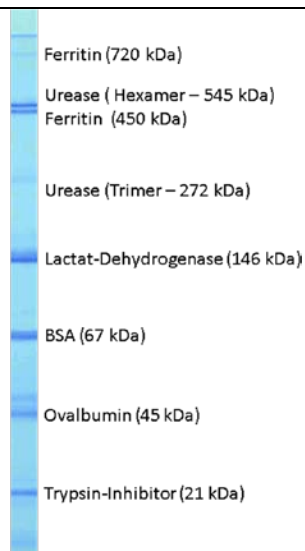

Storage

-20 °C, Storage in small aliquots is recommended to avoid repeated freeze/thaw cycles.

Operating instruction

Recommended loading volume:

- 5 – 10 µl/lane (gel size 10 x 8 cm, 0.75 oder 1 mm thickness)
- 10 – 15 µl/lane (gel size 30 x 20 cm, 1 oder 1.5 mm thickness)

Vers 0609